

Heelift – heel off-loading

3 versions for the best treatment

Effective off-loading for heel pressure ulcers, especially for bedridden patients

Heelift



- > The Original model – an ideal and economical device for heel off-loading

Heelift Glide



- > The premium version of the Heelift – providing an additional forefoot strap for more support and a gliding exterior

Heelift AFO



- > Heelift AFO – improved version combining all features of the Heelift Glide with a stabilizing plate for limited weight bearing


- > for the prevention and treatment of heel pressure ulcers, diabetic foot syndrome and decubitus
- > diabetes with neurological sensory deficit
- > prolonged immobility or coma
- > allows for patient-specific customization

2 types of interior foam available:
convoluted and smooth




DARCO

What makes the Heelift the right boot for your patient?

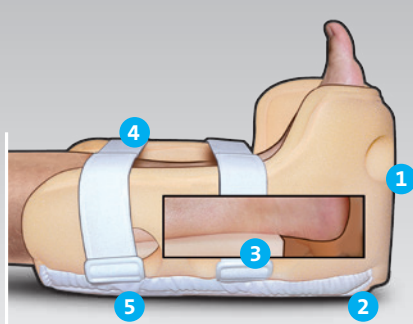


All Heelift products

- 1 **Ventilation Holes** provide additional air circulation to promote increased patient comfort.
- 2 **Heel Suspension Opening** allows for visual inspection of the heel.
- 3 **Affixed Elevation Pad** lifts the leg to suspend the heel, promoting the ideal environment for an ulcer to heal.
- 4 **Soft »easy-on, easy-off« Closure Straps** are a simple solution for a custom fit.
- 5 **Friction-free Tricot Backing** aids in patient mobility and helps ensure patient safety.

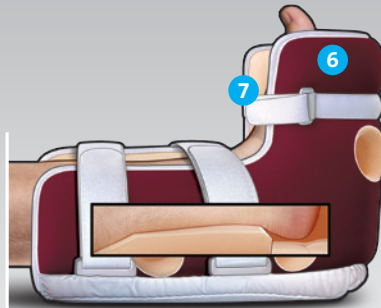

Heelift Glide

- 6 **Durable Exterior** allows the patient's leg to glide easily across bed sheets.
- 7 **Forefoot Strap** keeps the boot in place and helps support the forefoot.

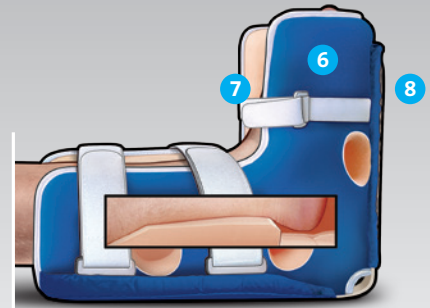

Heelift AFO

- 6 **Durable Exterior** allows the patient's leg to glide easily across bed sheets.
- 7 **Forefoot Strap** keeps the boot in place and helps support the forefoot.
- 8 **Semi-rigid brace** provides the best possible control of the ankle's position and motion.



Choice of two interiors

The **smooth** foam material was developed especially for patients with sensitive skin. This helps prevent the formation of oedema. This version is recommended for diabetic patients.

The **standard** version with convoluted foam material allows better air circulation. The attached nubs provide a pleasant cooling effect.

Item-No. Size (Calf circumference)

Heelift

Standard (convoluted foam)

DM-HLP Petite (15 – 25 cm)
DM-HL Normal (20 – 36 cm)

Smooth (smooth foam)

DM-HLSP Petite (15 – 25 cm)
DM-HLS Normal (20 – 36 cm)
DM-HLSB Bariatric (30 – 58 cm)
Colour: Beige

Heelift Glide

Standard (convoluted foam)

DM-HLPG Petite (15 – 25 cm)
DM-HLG Normal (20 – 36 cm)

Smooth (smooth foam)

DM-HLSPG Petite (15 – 25 cm)
DM-HLSG Normal (20 – 36 cm)
DM-HLSBG Bariatric (30 – 58 cm)
Colour: Red

Heelift AFO

Standard (convoluted foam)

DM-HLAFO Universal (20 – 36 cm)

Smooth (smooth foam)

DM-HLSAFO Universal (20 – 36 cm)
Colour: Blue

All Heelift products fit right and left. They can be washed safely at 50°C and can be autoclaved for sterility.

